

Alpaca anti-mouse IgG3, Fab-specific recombinant VHH, Alexa Fluor® 647

Product code: sms3AF647-1

Description	Monovalent, recombinant secondary single domain antibodies to mouse IgG3: Alpaca monoclonal Nanobody, Fab-specific, Alexa Fluor 647 conjugated
Product type	Secondary Nanobody
Format	Alpaca single domain antibodies, monovalent
Host	Alpaca-derived, recombinantly produced in bacteria
Target/ Specificity	This Nanobody mixture recognizes the Fab fragment of mouse IgG3.
Cross-reactivity	No cross-reactivity to rabbit, rat, sheep, goat, guinea pig, human, and macaque serum
Immunogen	Purified mouse IgG3
Clonality	Monoclonal
Clones	VHH0283
Conjugate	Site-directed conjugation to Alexa Fluor 647
Excitation/ Emission	Excitation max: 650 nm, Emission max: 665 nm
Synonyms	Alpaca single domain antibody, V _H H, Nanobody, binding domain of single domain antibody, Nano-antibody
Validation	Application validated for immunofluorescence and Western blotting Determination of cross-reactivity, subclass specificity, sequence, affinity, melting point, and degree of labeling (DOL)
Affinity (K_D)	VHH0283: K _D = 5.70 nM
DOL	3 fluorophores per Nanobody
Purity	Recombinantly expressed and purified
Form	Buffered aqueous solution
Storage buffer	10 mM HEPES pH 7.0, 500 mM NaCl, 5 mM EDTA, Preservative: 0.09% Sodium azide, Safety datasheet (SDS): Sodium azide SDS
Concentration	0.5 g/L
Size	10 µL; 100 µL
Storage instructions	Shipped at ambient temperature. Store at -20°C short term or -80°C long term. Aliquot upon delivery. Avoid freeze-thaw cycles. Stable for 6 months after arrival.
Applications	IF/ICC: recommended starting dilution 1:1000 Western blot: recommended starting dilution 1:1000 The optimal dilution depends on the application and should be determined by the user. A titration from range from 1:250 up to 1:2,000 is recommended. Note: Image acquisition time may have to be optimized.

Alpaca anti-mouse IgG3, Fab-specific recombinant VHH, Alexa Fluor® 647


Product code: sms3AF647-1

Tested applications

Immunofluorescence

Primary antibody: mouse anti-Lamin A/C antibody

Secondary antibody: alpaca anti-mouse IgG3, Fab-specific recombinant V_HH, Alexa Fluor 647 (sms3AF647-1) 1:1,000

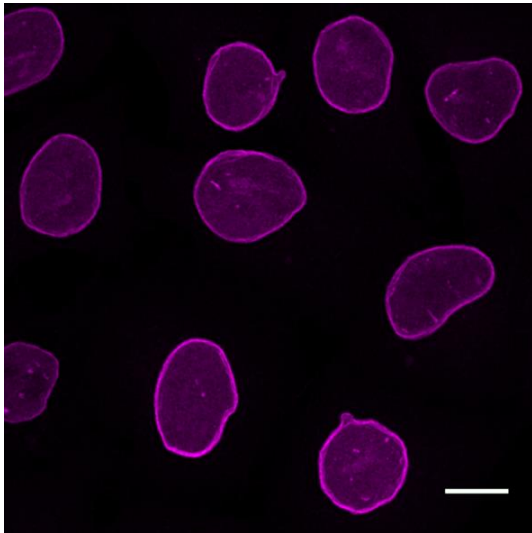

Immunostaining of Lamin A/C in HeLa cells with mouse anti-Lamin A/C antibody and alpaca anti-mouse IgG3, Fab-specific recombinant V_HH, Alexa Fluor 647 (magenta). Scale bar, 10 μ m.

Only for research applications, not for diagnostic or therapeutic use.

ChromoTek is a registered trademark of ChromoTek GmbH. Nanobody is a registered trademark of Ablynx, a Sanofi company. Alexa Fluor is a registered trademark of Life Technologies Corporation, a part of Thermo Fisher Scientific Inc.