

Alpaca anti-rabbit IgG, recombinant VHH, Alexa Fluor® 488

Product code: srbAF488-1

Description	Monovalent, recombinant secondary single domain antibodies to rabbit IgG: Mixture of 2 alpaca monoclonal Nanobodies, Fab- and Fc-specific, Alexa Fluor 488 conjugated
Product Type	Secondary Nanobody
Format	Alpaca single domain antibodies, monovalent
Host	Alpaca-derived, recombinantly produced in bacteria
Isotype	VHH domain of alpaca heavy chain antibody IgG2/3
Target/ Specificity	This Nanobody mixture recognizes Fab and Fc fragments of rabbit IgG.
Cross-Reactivity	No cross-reactivity to mouse, rat, sheep, goat, and guinea pig IgG/ serum Cross-reactivity to human and macaque IgG/ serum
Immunogen	Purified rabbit IgG
Clonality	Biclonal: mixture of 2 monoclonal Nanobodies
Clones	VHH0244, VHH0245
Conjugate	Site-directed conjugation to Alexa Fluor 488
Synonyms	Alpaca single domain antibody, VHH, Nanobody, binding domain of single domain antibody, Nano-antibody
Validation	Application validated for immunofluorescence and Western blotting Determination of cross-reactivity, sequence, affinity, melting point, and degree of labeling (DOL)
Affinity (K_D)	VHH0244: K _D = 0.18 nM, VHH0245: K _D = 1.2 nM
DOL	2 fluorophores per Nanobody
Purity	Recombinantly expressed and purified
Form	Buffered solution
Storage Buffer	10 mM HEPES pH 7.0, 500 mM NaCl, 5 mM EDTA, Preservative: 0.09% Sodium azide, Safety datasheet (SDS): Sodium azide SDS
Concentration	0.5 gram VHH per liter
Size	10 µL; 100 µL
Storage instructions	Shipped at ambient temperature. Store at -20°C short term or -80°C long term. Aliquot upon delivery. Avoid freeze-thaw cycles. Stable for 6 months after arrival.
Applications	IF/ICC: recommended starting dilution 1:1,000 Western blot: recommended starting dilution 1:1,000 The optimal dilution depends on the application and should be determined by the user. A titration from range from 1:250 up to 1:2,000 is recommended. Note: Image acquisition time may have to be optimized.


Alpaca anti-rabbit IgG, recombinant VHH, Alexa Fluor® 488

Product code: srbAF488-1

Tested applications

Immunofluorescence


Primary antibody: rabbit anti-Lamin B1 antibodies
Secondary antibody: alpaca anti-rabbit IgG, recombinant VHH, Alexa Fluor 488 (srbAF488-1) 1:1,000

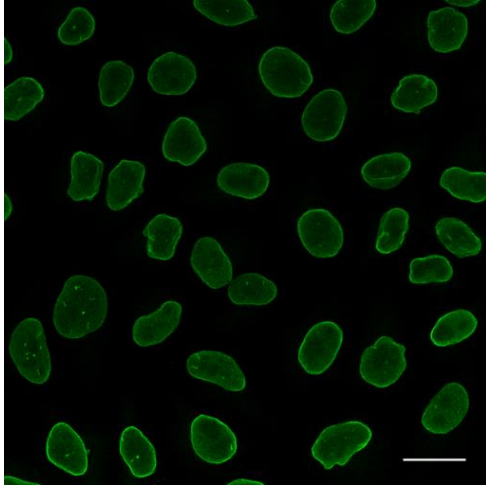

Immunostaining of nuclear lamina in HeLa cells with rabbit anti-Lamin B1 antibodies and alpaca anti-rabbit IgG, recombinant VHH, Alexa Fluor 488 (green). Scale bar, 20 μ m.

Western Blot

Primary antibody: rabbit anti-GFP PABG1 antibody (PABG1, ChromoTek) 1:1,000
Secondary antibody: alpaca anti-rabbit IgG, recombinant VHH, Alexa Fluor 488 (srbAF488-1) 1:1,000

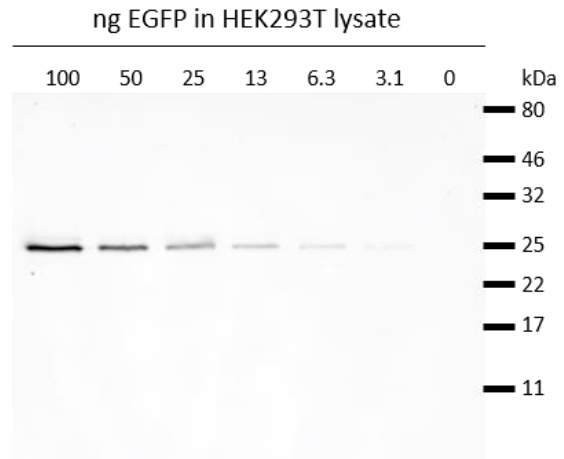

Western blot analysis of EGFP (EGFP-250, ChromoTek) added to HEK293T cell lysate. Detection with rabbit anti-GFP PABG1 antibody and alpaca anti-rabbit IgG, recombinant VHH, Alexa Fluor 488.

Only for research applications, not for diagnostic or therapeutic use.

ChromoTek is a registered trademark of ChromoTek GmbH. Nanobody is a registered trademark of Ablynx, a Sanofi company. Alexa Fluor is a registered trademark of Life Technologies Corporation, a part of Thermo Fisher Scientific Inc.