

Alpaca anti-mouse IgG1, Fc-specific recombinant VHH, Alexa Fluor® 568

Product code: sms1AF568-1

Description	Monovalent, recombinant secondary single domain antibodies to mouse IgG1: Mixture of 2 alpaca monoclonal Nanobodies, Fc-specific, Alexa Fluor 568 conjugated
Product type	Secondary Nanobody
Format	Alpaca single domain antibodies, monovalent
Host	Alpaca-derived, recombinantly produced in bacteria
Target/ Specificity	This Nanobody mixture recognizes the Fc fragment of mouse IgG1.
Cross-reactivity	No cross-reactivity to rabbit, rat, sheep, goat, guinea pig, human, and macaque serum
Immunogen	Purified mouse IgG1
Clonality	Biclonal: mixture of 2 monoclonal Nanobodies
Clones	VHH0302, VHH0305
Conjugate	Site-directed conjugation to Alexa Fluor 568
Excitation/ Emission	Excitation max: 578 nm, Emission max: 603 nm
Synonyms	Alpaca single domain antibody, V _H H, Nanobody, binding domain of single domain antibody, Nano-antibody
Validation	Application validated for immunofluorescence and Western blotting Determination of cross-reactivity, subclass specificity, sequence, affinity, melting point, and degree of labeling (DOL)
Affinity (K_D)	VHH0302: K _D = 0.1 nM, VHH0305: K _D = 0.6 nM
DOL	2 fluorophores per Nanobody
Purity	Recombinantly expressed and purified
Form	Buffered aqueous solution
Storage buffer	10 mM HEPES pH 7.0, 500 mM NaCl, 5 mM EDTA, Preservative: 0.09% Sodium azide, Safety datasheet (SDS): Sodium azide SDS
Concentration	0.5 g/L
Size	10 µL; 100 µL
Storage instructions	Shipped at ambient temperature. Store at -20°C/-4°F. Avoid freeze-thaw cycles. Aliquot upon arrival. Protect from light. Stable for 6 months.
Applications	IF/ICC: recommended starting dilution 1:500 (e.g. PBS supplemented with 4% BSA) Western blot: recommended starting dilution 1:1,000 (e.g. PBS supplemented with 0.075% Tween-20 and 5% skimmed milk) The optimal dilution depends on the application and should be determined by the user. A titration from range from 1:250 up to 1:2,000 is recommended. Note: Image acquisition time may have to be optimized.

Alpaca anti-mouse IgG1, Fc-specific recombinant V_HH, Alexa Fluor[®] 568


Product code: sms1AF568-1

Tested applications

Immunofluorescence

Primary antibody: mouse anti- α Tubulin antibody

Secondary antibody: alpaca anti-mouse IgG1, Fc-specific recombinant V_HH, Alexa Fluor 568 (sms1AF568-1) 1:500

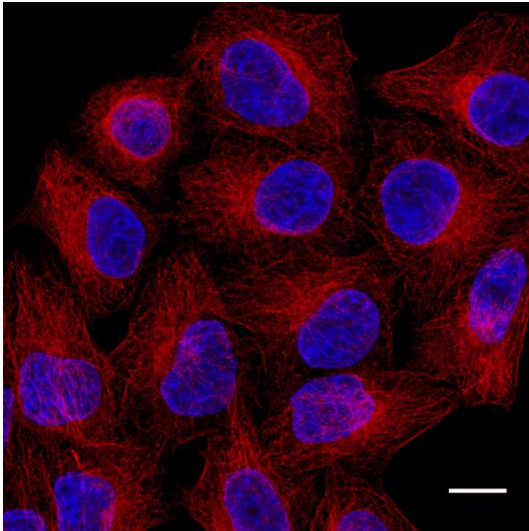


Immunostaining of Tubulin in HeLa cells with mouse anti- α Tubulin antibody and alpaca anti-mouse IgG1, Fc-specific recombinant V_HH, Alexa Fluor 568 (red). Nuclei were stained with DAPI (blue). Scale bar, 10 μ m.

Western Blot

Primary antibody: mouse anti-GFP antibody

Secondary antibody: alpaca anti-mouse IgG1, Fc-specific recombinant V_HH, Alexa Fluor 568 (sms1AF568-1) 1:1,000

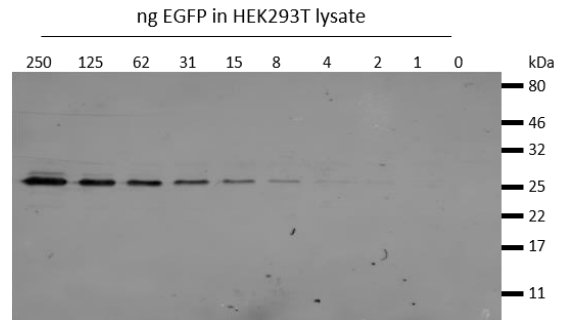

Western blot analysis of EGFP (EGFP-250, ChromoTek) added to HEK293T cell lysate. Detection with mouse anti-GFP antibody and alpaca anti-mouse IgG1, Fc-specific recombinant V_HH, Alexa Fluor 568.

Only for research applications, not for diagnostic or therapeutic use.

ChromoTek is a registered trademark of ChromoTek GmbH. Nanobody is a registered trademark of Ablynx, a Sanofi company. Alexa Fluor is a registered trademark of Life Technologies Corporation, a part of Thermo Fisher Scientific Inc.